



01-2010  
01-2020  
01-2011  
01-2021

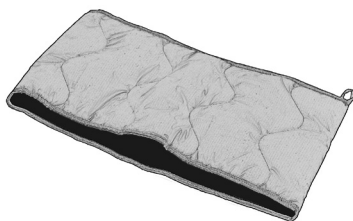

01-2015

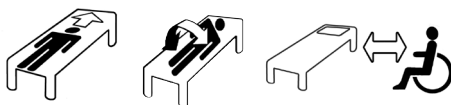

**SWL: 300 kg/660 lbs**

## Products

Order Ref.	Product	Description
01-2010	EasySlide	50 x 60 cm / 19,7 x 23,6 inch, Polyester
01-2020	EasySlide	50 x 60 cm / 19,7 x 23,6 inch, Microfibre
01-2011	EasySlide	50 x 70 cm / 19,7 x 27,2 inch, Polyester
01-2021	EasySlide	50 x 70 cm / 19,7 x 27,2 inch, Microfibre
01-2015	EasySlide	25 x 60 cm / 9,8 x 23,6 inch, Polyester

## Accessories

Order Ref.	Product	Description
01-2051	Protective cover	For EasySlide 01-2010, 01-2020, 01- 2030
01-2052	Protective cover	For EasySlide 01-2011, 01-2021
01-2055	Disposable cover	For EasySlide 01- 2010, 01- 2020, 01-2030

01-2030	SET	EasySlide 01-2020 + Protective cover 01- 2051
---------	-----	---



[info@amilakesouthern.co.uk](mailto:info@amilakesouthern.co.uk)

EasySlide is used for turning in bed, for transfers between bed and wheelchair/toilet, and for repositioning the patient higher up in bed. EasySlide is a tubular, padded mat with a friction-reducing material on the inside.



SystemRoMedic is a series of different models of transfer assistive devices. The choice of assistive devices depends on the transfer situation and the patient's functional capability and needs.

## Functional inspection

### Visual inspection

Check the condition and function of the assistive device regularly. Always inspect the product after laundering. Check to ensure that seams and material are free from damage. Check to ensure that fabric is not worn or faded. If there are signs of wear, the product must be discarded.

## Test friction

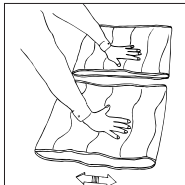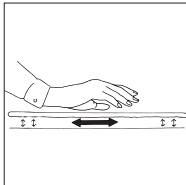

Test EasySlide's friction by rolling it and noting how easily it glides. If possible, compare with a new EasySlide.

## Always read the manual

Always read the manuals for all assistive devices used during a transfer. Keep the manual where it is accessible to users of the product. Do not leave the patient unattended during a transfer situation.

## Placement of EasySlide

Place the EasySlide under one buttock or under one shoulder. For transfers, all of the patient's body weight should be shifted over onto the EasySlide.


To remove, pull on the lower end of the EasySlide.

## Using EasySlide



1. EasySlide can be used to advantage for transfers to/from a wheelchair. Complement, preferably, with an EasyGlide.

2. EasySlide can be used under the shoulders as a friction-reducing assistive device to facilitate repositioning higher up in bed.

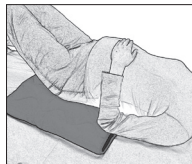

3. EasySlide can be used under the hips to facilitate lateral transfer. Active patients can use EasySlide themselves, but it can also be used by less active patients with the aid of a caregiver.



The gap between the underlying surfaces (bed/gurney/show trolley) must not exceed 5 cm.

NOTE!

Always remember to lock wheels.



To prevent the risk of the patient sliding/falling, never leave the EasySlide under the patient after the completed transfer, and never position the patient too close to the edge of the bed.

## Material:

Polyester, microfibre, polyamide

## Care of the product:

Read the product label.



Do not use rinsing agent.

To prolong product life, avoid tumble-drying.

**CE** Medical Device Class I

The product meets the requirements of the Medical Devices Directive 93/42/EEC.

# SystemRoMedic™

Easy-to-use assistive devices for

Easy transfers and to making life and work easier for both patients and personnel in the care sector.

Experience, innovation and training are the basis for SystemRoMedic, a total solution for every imaginable transfer situation.

**Transfer:** products for moving patients between locations.

**Positioning:** products for repositioning in the same location.

**Support:** products that provide mobility support.

**Lifting:** products adapted for lifting.

The philosophy behind SystemRoMedic is to prevent occupational injuries while improving the patient's sense of independence and dignity. Through a combination of training and a complete range of transfer-assistive devices, SystemRoMedic offers the means to improving both the work environment and the quality of care while enabling significant cost savings.

Their mission, is to help people, has always been, and will continue to be, their driving force of innovation. They love easy transfers.



[info@amilakesouthern.co.uk](mailto:info@amilakesouthern.co.uk)